

## VETPUMP® 2

Precision Suction &  
Irrigation System



# Cleaning ears has never been so easy

The VETPUMP® 2 from KARL STORZ puts the safest and most effective flushing and suction at your fingertips!

More powerful, more precise and safer than other units, the VETPUMP® 2 is the ultimate device for lavage and suction during otoscopy, gastroscopy and more.



## Features:

- Powerful suction and irrigation pump for both rigid and flexible KARL STORZ endoscopes
- There is no exchange of air between bottles, eliminating any potential contamination of clean irrigation fluids
- Micro-processor controlled unit is quiet and responsive; unit is completely silent when trumpet valves are not activated
- Irrigation pressure and suction vacuum can be regulated separately
- User-selected pressures are precise, powerful and reproducible
- Versatile hand piece attaches to catheters, irrigation tubes, or directly to endoscope channel
- Bottles, tubing and handle are fully autoclavable

Check out our schedule of upcoming hands-on training courses at <http://go.karlstorz.com/eventsVET>





69321601

69321601

**VETPUMP® 2**, power supply 100/120/230/240 VAC, 50/60 Hz including:



**Suction Bottle**, 1.5 l, sterilizable  
**Bottle Cap**, for use with suction bottle 1.5 l, sterilizable




**Irrigation Bottle**, 1 l, sterilizable  
**Bottle Cap**, for use with irrigation bottle

**Bottle Holder**, for suction bottle  
**Bottle Holder**, for irrigation bottle  
**Bacterial Filter**

**Mains Cord**  
**Mains Cord**, US-version

Pressure Regulated (mmHg)	Flow Regulated (l/min)	Suction Pressure (mmHg)	Power Supply	Dimensions (w x h x d) mm	Weight (kg)	Certified to
max. 600	max. 3.5	regulated: max. 580	100/120/230/240 VAC, 50/60 Hz	305 x 110 x 260	3.8	EN 61010-1

## Accessories for Otoscopy:

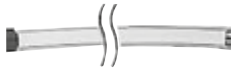

69320033 **Silicone Tubing Set**, irrigation-suction from unit to bottles, **autoclavable**



69811 **Instrument Tubing Set**, autoclavable, to connect VETPUMP® 2 bottles to handle



69810 **Handle**, for irrigation and suction, autoclavable, for use with tubing set 69811, in conjunction with 5 mm suction and irrigation tubes or catheters any size with LUER Lock connection



## Accessories for Abdominal/Thoracic Lavage:



69320033 **Silicone Tubing Set**, irrigation-suction from unit to bottles, **autoclavable**



69811 **Instrument Tubing Set**, autoclavable, to connect VETPUMP® 2 bottles to handle



37112A **Handle**, straight, with clamping valve, for suction and irrigation, **autoclavable**, for use with 5 and 10 mm suction and irrigation tubes



69811LAP **Adaptor**, for use with suction and irrigation Handle 37112A, plus Instrument Tubing Set 69811, package of 5



031121-10 **Gas filter**, for single use, sterile, package of 10, with male LUER-Lock connector for attachment to irrigation connector on unit



## Accessories for Gastrointestinal Endoscopy:



20300042 **Silicone Tubing Set**, for suction only, sterilizable




610VE **Silicone Tubing Set**, for insufflation, attaches to light source adaptors 487VE or 69487KVE




# VETPUMP® 2

## System Components for use in Otoscopy

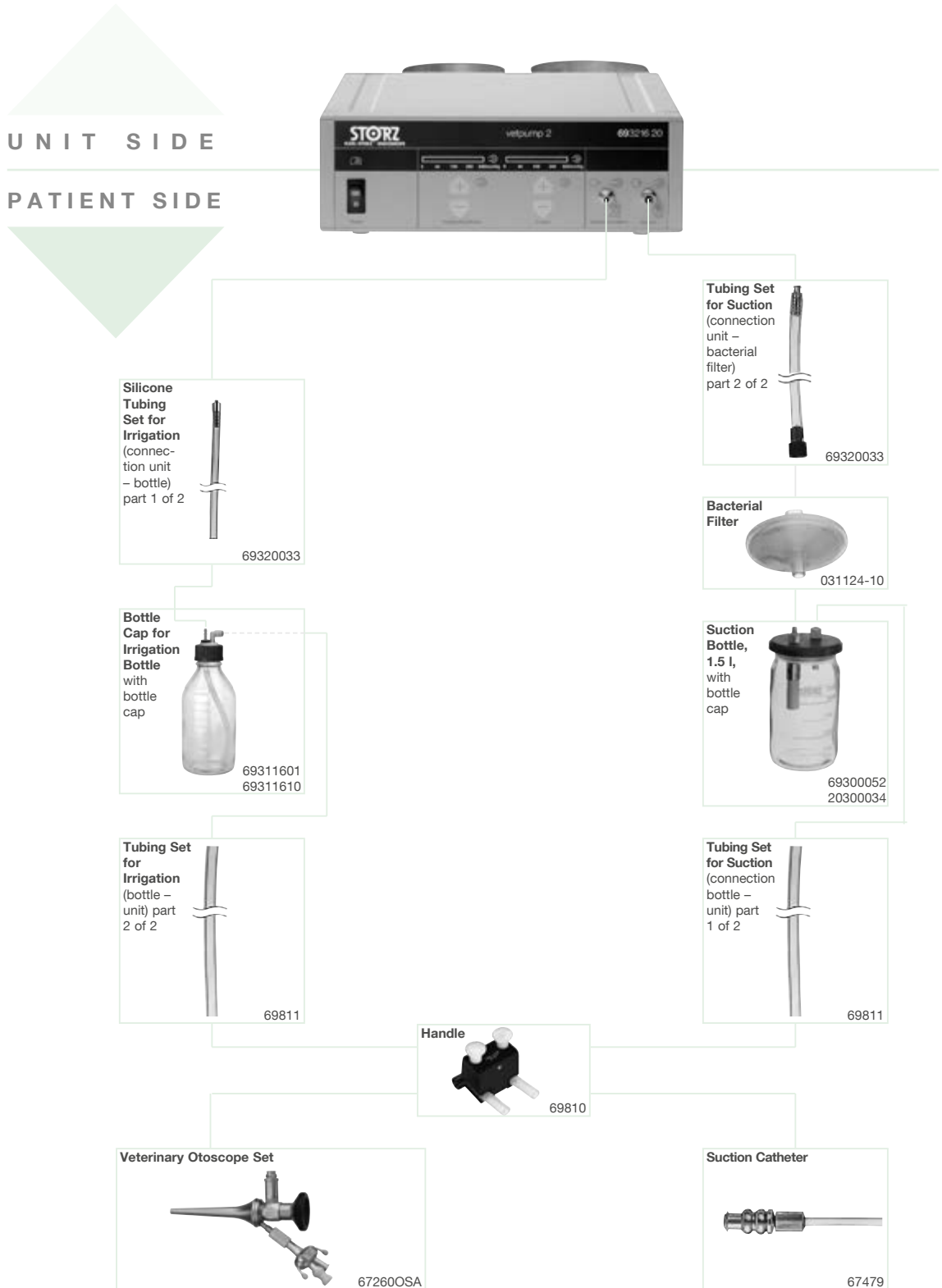

# VETPUMP® 2


## System Components for use in Laparoscopy / Thoracoscopy

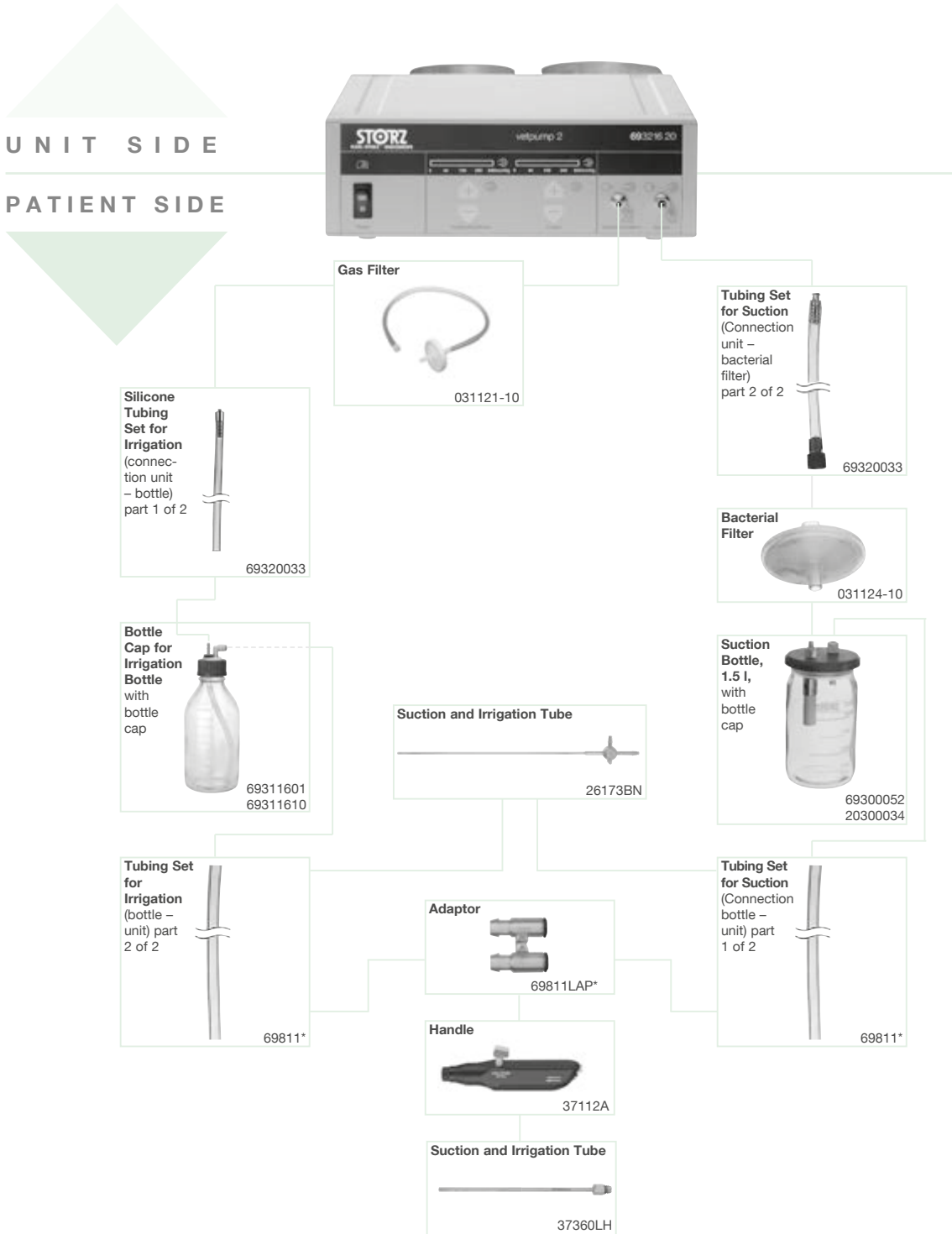

\* Remove black plastic fitting from Tubing Set 69811 and replace with Adaptor 69811LAP.

# VETPUMP® 2

## System Components for use in Gastrointestinal Endoscopy

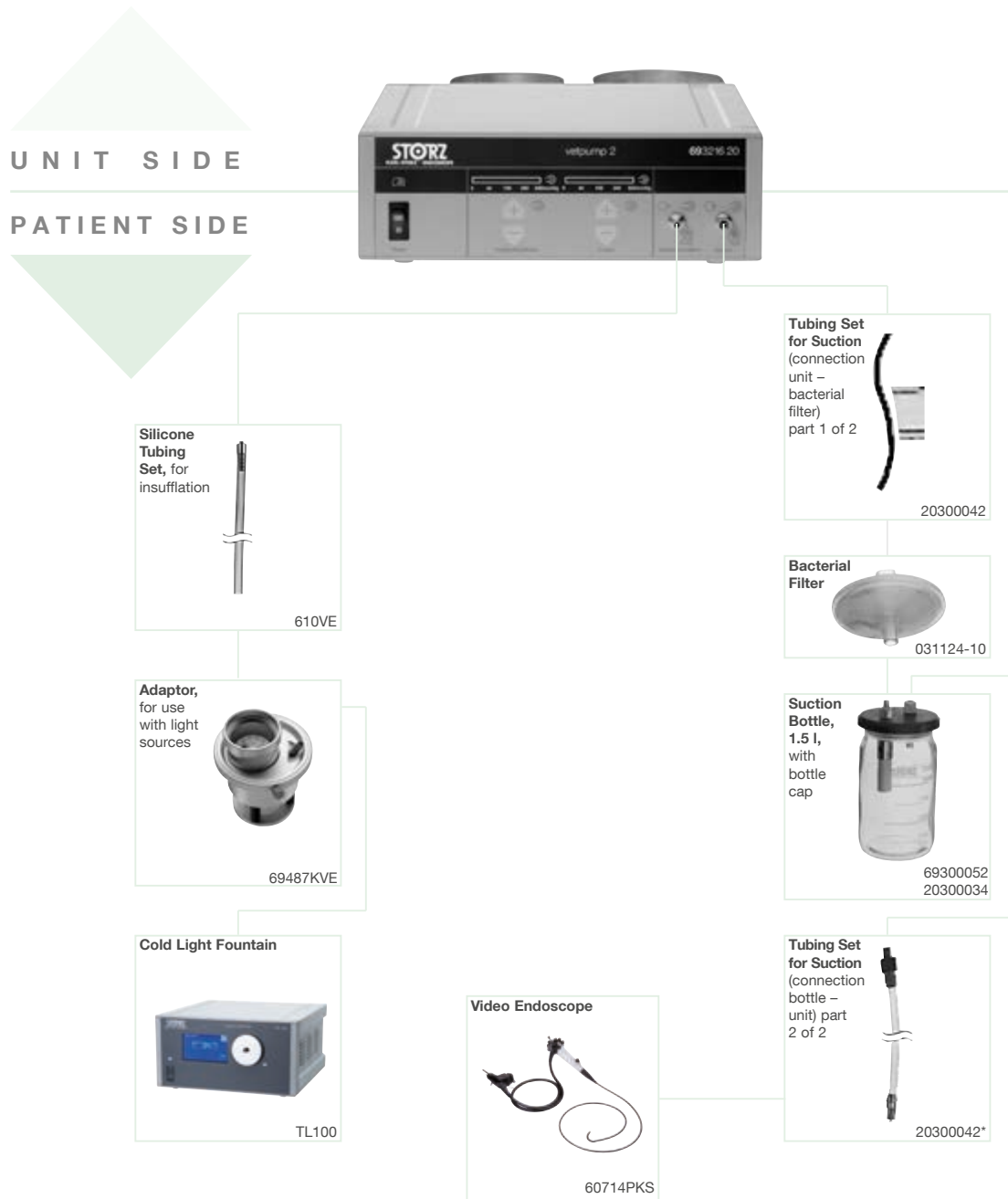

\* Remove male LUER Lock fitting from end of Tubing 20300042 and connect free end to suction connector of endoscope.

More than  
**75**  
Years

*Shaping the Future  
of Endoscopy with you*

**STORZ**  
KARL STORZ—ENDOSKOPE

THE DIAMOND STANDARD

KARL STORZ SE & Co. KG  
Dr.-Karl-Storz-Straße 34, 78532 Tuttlingen/Germany  
Postbox 230, 78503 Tuttlingen/Germany  
Phone: +49 7461 708-0  
Fax: +49 7461 708-105  
E-Mail: [vet@karlstorz.com](mailto:vet@karlstorz.com)  
[www.karlstorz.com](http://www.karlstorz.com)